HOMECARE WEIGHING SOLUTION™

Body Fat Analyzer

# GBF-830

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## CAS

[www.cas.co.kr](http://www.cas.co.kr)

OWNER'S MANUAL

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**Thank you for your purchase.**

For the proper using & maintaining the scale please read this instruction completely before operating.

**TABLE OF CONTENT**

LCD FULL DISPLAY ..... 11

FEATURES ..... 11

FUNCTION KEYS ..... 11

TO CHANGE BATTERIES ..... 12


TO INPUT PERSONAL INFORMATION ..... 12

TO WEIGH AND MEASURE ..... 13

BODY FAT, WATER AND MUSCLE RATE RANGE ..... 15


## LCD FULL DISPLAY
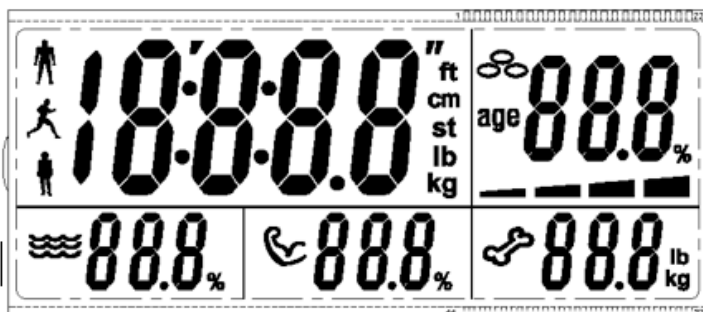
## FEATURES

- Capacity: 180 kg
- Graduation: 0.1 kg
- Auto on or press Unit to switch on
- 3 function key: SET, ▲, ▼.
- It may measure body fat, water, muscle and bone
- Identify automatically function: can store 8 users' parameter and the last measure data
- Auto on and Auto off function
- Over load indication "- - - -"
- Low battery indication: "Lo"
- Power: 4xAAA Battery


## FUNCTION KEYS
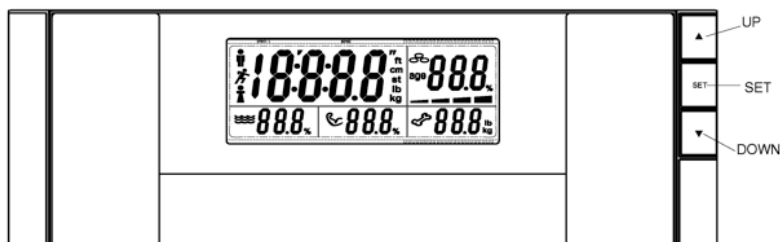
## TO CHANGE BATTERIES

- Ensure +/- terminals are correct.
- Use proper tool to help release the battery from the compartment if necessary.
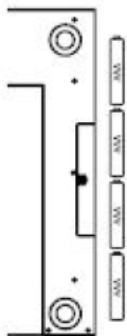
## TO INPUT PERSONAL INFORMATION

- Set up personal data(gender/age/stature) before using, it can set up 8 users' data
1. Press "set" to enter the mode of setup, the flashing Pn (n can be 1~8)means the memory of the data.  
Press ▲&▼ to change number, gender, stature, age, then press "SET" to confirm. The scale will back to the mode of weighing after finishing setting up


2. Set up the gender: press UP, DOWN (can select female, male, athlete mode), then press SET to confirm

3. Set up the stature: Press UP, DOWN, and then press SET to confirm

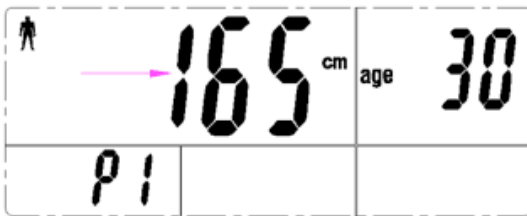4. Set up age: Press UP, DOWN, and then press SET to confirm

## TO WEIGH AND MEASURE

1. The weighing after the personal data has been set in the scale

1-1 After setting up of the data, when LCD displays 0.0, you may weigh or stand on the scale platform directly, the weight value will be shown after a few seconds.

1-2 After the weight is stable, the flashing "0" means that the scale is measuring the fat, the "0" signal moves from left to right for two times, and then shows the data

1-3 After the measuring, the LCD will display the weight, fat, water, muscle and bone data, after 20+/-4S the scale will be automatically turned off.

2. The weighing without being set up personal data in the scale


2-1 After the weight is stable, if the scale can't identify the user, the weight will be shown on the scale for 6+/-1.5s and then automatically off.

2-2 If the weight is stable and can identify the user, the scale can measure the body fat

2-3 After "0" stops flashing, the scale can identify the user by weight, and press SET to select the user.

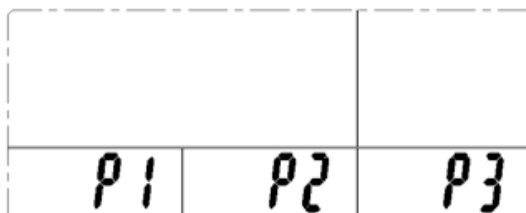2-4 After the user is selected, the LCD displays weight, fat, water, muscle, bone data, and after 20+/-4S the scale is automatically off.

### BODY FAT, WATER AND MUSCLE RATE RANGE

<1> The body fat ratio (%): 3%-50%

Age	Female				Male			
	Underfat	Healthy	Overfat	Obese	Underfat	Healthy	Overfat	Obese
<39	<21	21.1-33	33.1-39	>39	<8	8.1-20	20.1-25	>25
40-59	<23	23.1-34	34.1-40	>40	<11	11.1-22	22.1-28	>28
>59	<24	24.1-36	36.1-42	>42	<13	13.1-25	25.1-30	>30

<2> The body water ratio(%): 25%-75%

Gender	Female			Male		
	Low	Normal	High	Low	Normal	High
Age						
18-30	<47	47.1-57	>57	<53	53.1-67	>67
31-60	<42	42.1-52	>52	<47	47.1-61	>61
61-80	<37	37.1-47	>47	<42	42.1-56	>56

<3> The body muscle ratio(%): 25%-75%

Gender	Female	Male
Age	Healthy	Healthy
18~100	>34	>40

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