Thank you for your purchase.

Thank you for choosing this electronic personal scale of our series products. To ensure proper operation and to maintain the lifespan and the utmost accuracy of the scale, please read this manual carefully before use.

Key Features


- Capacity: 6kg-150kg
- Division: 100g
- Auto on & Auto off function
- Low battery and over load indication
- Free switch within three different units: kg/lb/st
- Powered by 4*AAA are packed in the box

Weighing Steps

1. If it is the first time to use, please take out 4xAAA batteries from insulation sheet, open the battery cover at the back of the scale. Put the 4xAAA batteries into the scale, and then close the cover.
2. Make sure to place scale on the flat floor.
3. Stand on the center of the weighing platform and keep still, LCD will show the result flash.
4. Come down from the scale when finish.
5. Scale will auto power off within 15 seconds after the weighing result is stable.

Weighing Steps

- LCD display means overload.

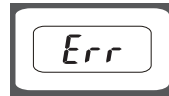- LO means power not enough.

**Attention**

1. Capacity: 6kg-150kg, d=100g, Please weigh within this scope, otherwise, the sensor will be damaged.
2. Do not put anything on the glass after weighing, or will short the life of battery.
3. Only for body weight measure, do not use elsewhere.
4. Do not use chemical detergent for cleaning, soft silk will be good.
5. Please stand on the center of the scale, keep you feet and the platform surface dry, avoid falling down and slip.
6. Do not put the scale on rough floor or carpet weighing, otherwise, the measure result will not precise.
7. Keeping the scale from wet, hot and cold place.
8. Do not jump or rough stomp on scale, not attempt to Disassembly.
9. If it can't work, please check whether battery touched electrodes or replace battery, or contact local agents for help.

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CAS
www.globalcas.com
OWNER'S MANUAL

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