HEALTHCARE SOLUTION™

Digital Personal Scale

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HE-53

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Thank you for your purchase.

For the proper using & maintaining the scale please read this instruction completely before operating.


TABLE OF CONTENT

SPECIFICATIONS	8
HOW TO USE	8
FUNCTION	
- ERROR INDICATION	8
- UNIT EXCHANGE(kg / lb)	9
- INSERT THE BATTERY	9
ADVICE FOR USE AND CARE	10

SPECIFICATIONS

- Equipped with a high precision strain gauge sensor system
- Capacity: 150kg (330lb)
- Division:0.05kg/0.1lb
- Safety and environmental Fire-proof platform
- Mini effective weight: 8kg
- Auto zero resetting
- Auto on / off
- Over-load indication
- Power: 1 x 3V CR2032 lithium battery

HOW TO USE

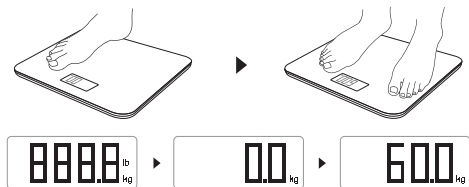- Put the scale on a hard & flat surface.
- Step on the scale steadily, till the weight displays on the screen. The scale will power off automatically after 10 seconds no operation.

<NOTICE> Please retry the scale when the weight is a bit different.

FUNCTION

ERROR INDICATION

- **Over-load Indication**

When the scale is over-loaded (more than 150kg), the screen will appear "Err".

• Low battery indication

When the screen indicates “Lo”, it will turn off automatically after 2 seconds. Open the battery cover and replace with new batteries in correct polarity.


UNIT EXCHANGE(kg/lb)

Under weight mode, press the unit key to choose the unit, the key is on the back of the scale.

Note: Please set the unit before weighing. It could not change the unit while weighing.

INSERT THE BATTERY

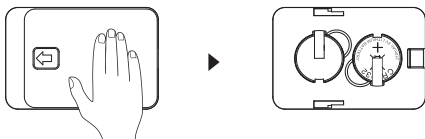- Pull the battery cover as per the arrow direction.
- Install 1 x 3V CR2032 lithium battery as shown.
- After installing the battery, put back the battery cover and push the cover backward.
- Please be aware that the incorrect polarity will cause the permanent impact on the electronic parts.

ADVICE FOR USE AND CARE

- For accurate results, place scale on a flat surface.
- Static pressure is effective, do not drop it or strike the scale, otherwise the scale is easy to damage.
- In weighing, please keep balance. Do not let the body in wobbling.
- Please use a damp cloth to clean the scale. Do not allow water to get inside scale.
- Do not use chemicals to clean the scale.
- To ensure the life of the scale, do not put the scale under wet & hot condition.
- Do not use for legal trade.
- Do not overload, otherwise the sensor will be damaged.
- Do not put heavy things on the scale when not use.
- Do not store the scale in an upright position.
- Due to the continuous improvements, this product is subject to change without notice.
- Please do not jump or stamp on the scale and avoid dropping it, as the scale is a highly precise electronic instrument
- Do not use it on the wet and slippery floor, do not stand on at the edge of scale to avoid any falling down.

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